

Alpha SureFire® Ultra™ Multiplex

Human and Mouse p-NF-κB p65 (Ser536) + Total NF-κB p65 Detection Kit

Product number: MPSU-PTNFKB-K500, MPSU-PTNFKB-K10K,
MPSU-PTNFKB-K50K, MPSU-PTNFKB-K-HV



Kit specificity:


This assay kit contains antibodies which recognize the phospho-Ser536 epitope and a distal epitope on NF-κB, plus antibodies to detect Total NF-κB p65. The protein detected in this kit corresponds to UniProt ID Q04206. NF-κB p65 is also known as Transcription factor p65, RELA and NFKB3. These antibodies recognize NF-κB p65 of human and mouse origin. Other species should be tested on a case-by-case basis.

Control lysate information:

Positive Control Lysate: Prepared from HeLa cells, cultured to confluence in T175 flasks in 10% FBS containing medium, treated with 50 ng/mL TNFα and 20 ng/mL Calyculin A for 10 minutes and lysed with 10 mL of Lysis Buffer.

Representative data

Data obtained from a 2-plate, 2-incubation protocol. HeLa cells were seeded at 20K cells/well in a 96 well plate and incubated overnight. Cells were serum starved for 2 hours, then treated with TNFα at the indicated concentrations for 5 minutes. Cells were lysed with Lysis Buffer and assayed for Phospho (536) and Total NF-κB p65 in the same wells of an assay plate using the SureFire Ultra Multiplex kit. Equivalent to approximately 1,000 cells/datapoint.

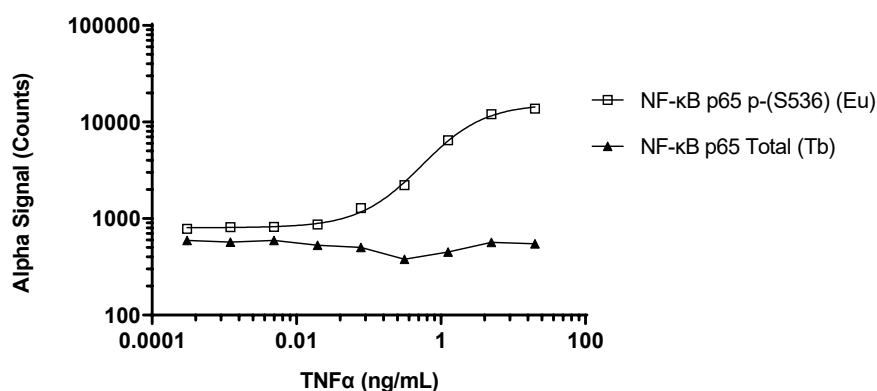

The information provided in this document is for reference purposes only and may not be all-inclusive. Revvity, Inc., its subsidiaries, and/or affiliates (collectively, "Revvity") do not assume liability for the accuracy or completeness of the information contained herein. Users should exercise caution when handling materials as they may present unknown hazards. Revvity shall not be liable for any damages or losses resulting from handling or contact with the product, as Revvity cannot control actual methods, volumes, or conditions of use. Users are responsible for ensuring the product's suitability for their specific application. REVVITY EXPRESSLY DISCLAIMS ALL WARRANTIES, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, REGARDLESS OF WHETHER ORAL OR WRITTEN, EXPRESS OR IMPLIED, ALLEGEDLY ARISING FROM ANY USAGE OF ANY TRADE OR ANY COURSE OF DEALING, IN CONNECTION WITH THE USE OF INFORMATION CONTAINED HEREIN OR THE PRODUCT ITSELF"