

AlphaLISA® SureFire® Ultra™

Human and Mouse p-c-Jun (Ser73) Detection Kit

Product number: ALSU-PCJUN-B500, ALSU-PCJUN-B10K,
 ALSU-PCJUN-B50K, ALSU-PCJUN-B-HV



Kit specificity:


This assay kit contains antibodies which recognize the phospho-Ser73 epitope and a distal epitope on c-Jun. The protein detected by this kit corresponds to UniProt ID P05412. c-Jun is also known as transcription factor c-Jun, AP1 and p39. These antibodies recognize c-Jun of human and mouse origin. Other species should be tested on a case-by-case basis.

Control lysate information:

Positive Control Lysate: Prepared from NIH/3T3 cells, cultured to confluence in T175 flasks in 10% FBS containing medium, treated with 500 ng/mL Anisomycin for 30 minutes and lysed with 7 mL of Lysis Buffer.

Representative data

Data obtained with a 2-plate, 2-incubation protocol. NIH/3T3 cells were seeded at 40K cells/well in a 96 well plate and grown overnight. Cells were treated with Anisomycin at the indicated concentrations for 30 minutes. Cells were lysed with Lysis Buffer and assayed separately for various targets using respective *SureFire Ultra* kits. Equivalent to approximately 2,000 cells/datapoint.

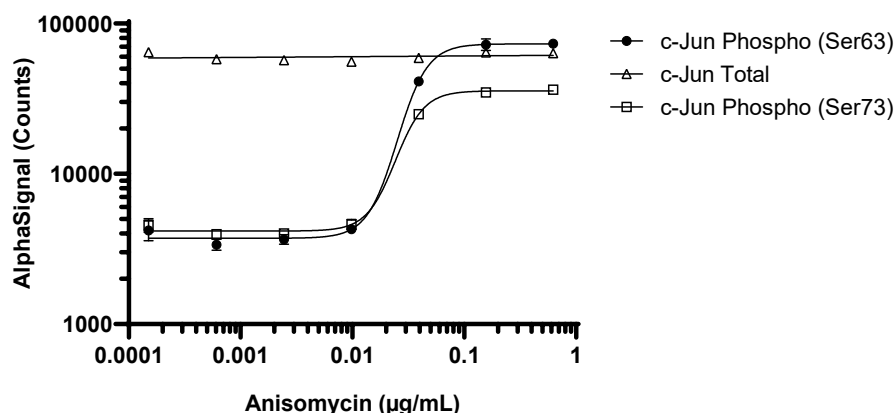

The information provided in this document is for reference purposes only and may not be all-inclusive. Revvity, Inc., its subsidiaries, and/or affiliates (collectively, "Revvity") do not assume liability for the accuracy or completeness of the information contained herein. Users should exercise caution when handling materials as they may present unknown hazards. Revvity shall not be liable for any damages or losses resulting from handling or contact with the product, as Revvity cannot control actual methods, volumes, or conditions of use. Users are responsible for ensuring the product's suitability for their specific application. REVVITY EXPRESSLY DISCLAIMS ALL WARRANTIES, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, REGARDLESS OF WHETHER ORAL OR WRITTEN, EXPRESS OR IMPLIED, ALLEGEDLY ARISING FROM ANY USAGE OF ANY TRADE OR ANY COURSE OF DEALING, IN CONNECTION WITH THE USE OF INFORMATION CONTAINED HEREIN OR THE PRODUCT ITSELF"