

AlphaLISA® SureFire® Ultra™

Mouse p-PERK (Thr980) Detection Kit

Product number: ALSU-PPERK-A500, ALSU-PPERK-A10K,
ALSU-PPERK-A50K, ALSU-PPERK-A-HV



Kit specificity:


This assay kit contains antibodies which recognize the phospho-Thr980 epitope and a distal epitope on PERK. The protein detected by this kit corresponds to UniProt ID Q9Z2B5. PERK is also known as EIF2AK3 and eukaryotic translation initiation factor 2 alpha kinase 3. These antibodies recognize PERK of mouse origin. Other species should be tested on a case-by-case basis.

Control lysate information:

Positive Control Lysate: Prepared from L-929 cells, cultured to confluence in T175 flasks in 10% FBS containing media, then treated with 500 nM Thapsigargin for 20 minutes and lysed with 4 mL of Lysis Buffer.

Representative data

Data obtained with a 2-plate, 2-incubation protocol. L929 cells were plated at 40K cells/well in a 96 well plate and incubated overnight. Cells were treated with Thapsigargin at the indicated concentrations for 20 minutes. Cells were lysed with Lysis Buffer and assayed separately for Phospho (Thr980) and Total PERK, using respective *SureFire Ultra* kits. Equivalent to approximately 8,000 cells/datapoint.

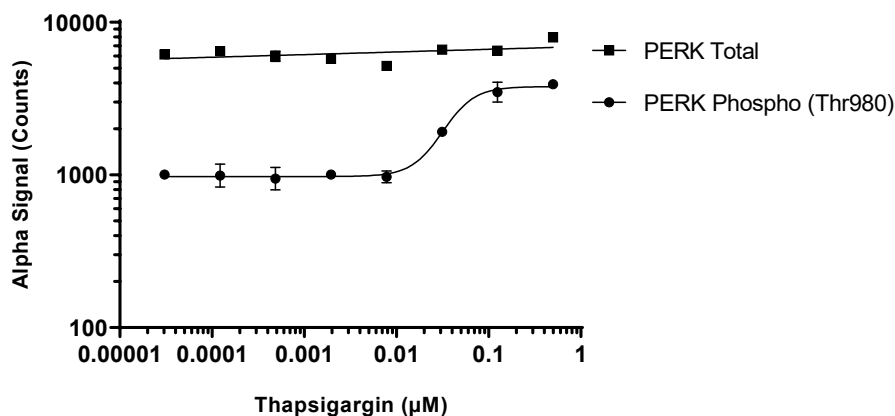

The information provided in this document is for reference purposes only and may not be all-inclusive. Revvity, Inc., its subsidiaries, and/or affiliates (collectively, "Revvity") do not assume liability for the accuracy or completeness of the information contained herein. Users should exercise caution when handling materials as they may present unknown hazards. Revvity shall not be liable for any damages or losses resulting from handling or contact with the product, as Revvity cannot control actual methods, volumes, or conditions of use. Users are responsible for ensuring the product's suitability for their specific application. REVVITY EXPRESSLY DISCLAIMS ALL WARRANTIES, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, REGARDLESS OF WHETHER ORAL OR WRITTEN, EXPRESS OR IMPLIED, ALLEGEDLY ARISING FROM ANY USAGE OF ANY TRADE OR ANY COURSE OF DEALING, IN CONNECTION WITH THE USE OF INFORMATION CONTAINED HEREIN OR THE PRODUCT ITSELF"