

AlphaLISA® SureFire® Ultra™

Human and Mouse p-STAT2 (Tyr690) Detection Kit

Product number: ALSU-PST2-A500, ALSU-PST2-A10K,
ALSU-PST2-A50K, ALSU-PST2-A-HV



Kit specificity:

This assay kit contains antibodies which recognize the phospho-Tyr690 epitope and a distal epitope on STAT2. The protein detected by this kit corresponds to UniProt ID P52630. STAT2 is also known as Signal transducer and activator of transcription 2. These antibodies recognize STAT2 of human and mouse origin. Other species should be tested on a case-by-case basis.


Control lysate information:

Positive Control Lysate: Prepared from A431 cells, cultured to confluence in T175 flasks in 10% FBS containing medium, then treated with 10 ng/mL recombinant human IFN β for 30 minutes and lysed with 4 mL of Lysis Buffer.


Representative data:

Data obtained with a 2-plate, 2-incubation protocol. RPMI 8226 cells were harvested, seeded in a 96-well plate at 200K cells/well and treated with 10 ng/mL IFN α at the indicated timepoints. Cells were washed, lysed with Lysis Buffer and assayed separately for Phospho (Tyr690) and Total STAT2 using respective *SureFire Ultra* kits. Equivalent to approximately 20,000 cells/datapoint.


RPMI 8226 cells treated with IFN α

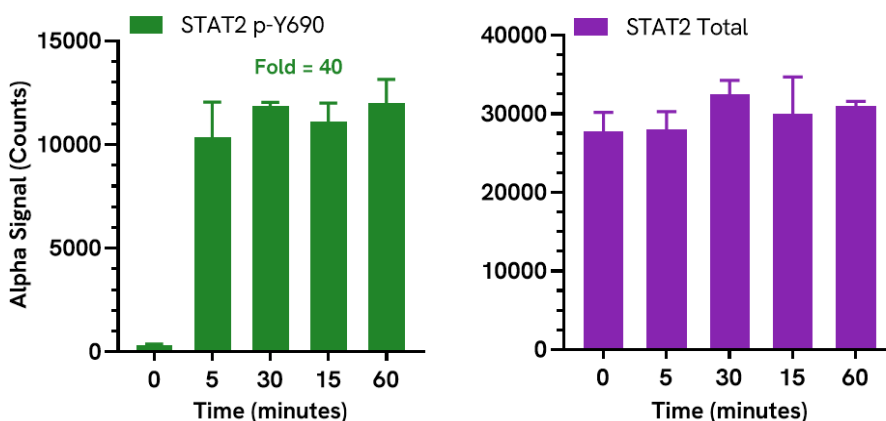

Data obtained with a 2-plate, 2-incubation protocol. A431 cells were seeded at 60K cells/well in a 96-well plate and incubated overnight. Cells were treated with IFN α at the indicated concentrations for 30 minutes. Cells were lysed with Lysis Buffer and assayed separately for Phospho (Tyr690) and Total STAT2 using respective *SureFire Ultra* kits. Equivalent to approximately 6,000 cells/datapoint.

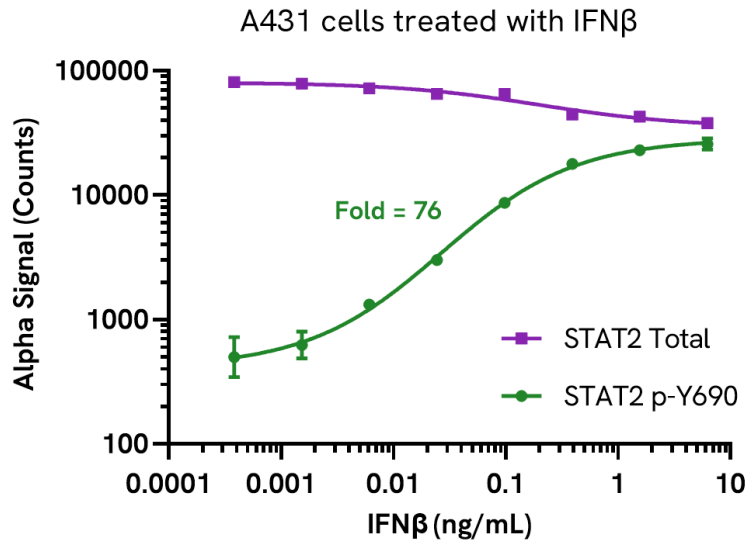

Data obtained with a 2-plate, 2-incubation protocol. THP-1 cells were seeded at 200K cells/well in a 96-well plate with HBSS containing 10 ng/mL IFN β for 15 minutes. Cells were then treated with the JAK inhibitor, Ruxolitinib at the indicated concentrations for 1 hour. Cells were lysed with Lysis Buffer and assayed separately for Phospho (Tyr690) and Total STAT2 using respective *SureFire Ultra* kits. Equivalent to approximately 20,000 cells/datapoint.

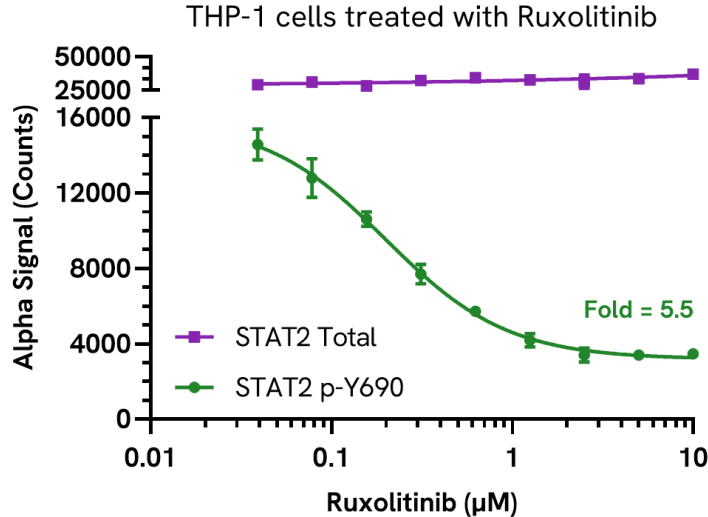

The information provided in this document is for reference purposes only and may not be all-inclusive. Revvity, Inc., its subsidiaries, and/or affiliates (collectively, "Revvity") do not assume liability for the accuracy or completeness of the information contained herein. Users should exercise caution when handling materials as they may present unknown hazards. Revvity shall not be liable for any damages or losses resulting from handling or contact with the product, as Revvity cannot control actual methods, volumes, or conditions of use. Users are responsible for ensuring the product's suitability for their specific application. REVVITY EXPRESSLY DISCLAIMS ALL WARRANTIES, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, REGARDLESS OF WHETHER ORAL OR WRITTEN, EXPRESS OR IMPLIED, ALLEGEDLY ARISING FROM ANY USAGE OF ANY TRADE OR ANY COURSE OF DEALING, IN CONNECTION WITH THE USE OF INFORMATION CONTAINED HEREIN OR THE PRODUCT ITSELF"

