

AlphaLISA™ SureFire® Biotin Free

Human CDK4 Total Detection Kit

Product number: ASBF-TCDK4-A500, ASBF-TCDK4-A10K,
ASBF-TCDK4-A50K, ASBF-TCDK4-A-HV



Kit specificity:


This assay kit contains antibodies that recognize distinct epitopes on CDK4. The protein detected by this kit corresponds to UniProt ID P11802. CDK4 is also known as Cyclin-dependent kinase 4. These antibodies recognize CDK4 of human origin. Other species should be tested on a case-by-case basis.

Control lysate information:


Positive Control Lysate: Prepared from HeLa cells, cultured to confluence in T175 flasks in 10% FBS containing media and lysed with 5 mL of Lysis Buffer.

Representative data:

Data obtained with a 2-plate, 2-incubation protocol. MCF7 cells were seeded at 40K cells/well in a 96-well plate and incubated overnight. Cells were treated with 5μM LY294002 (PI3K inhibitor) for 3 hours and then with increasing concentrations of Neuregulin-1 (NRG-1) for 24 hours in serum free media. Cells were lysed with Lysis Buffer and assayed separately for Total CDK4 and ERK using respective SureFire Biotin Free kits. Equivalent to approximately 2,000 cells/datapoint.

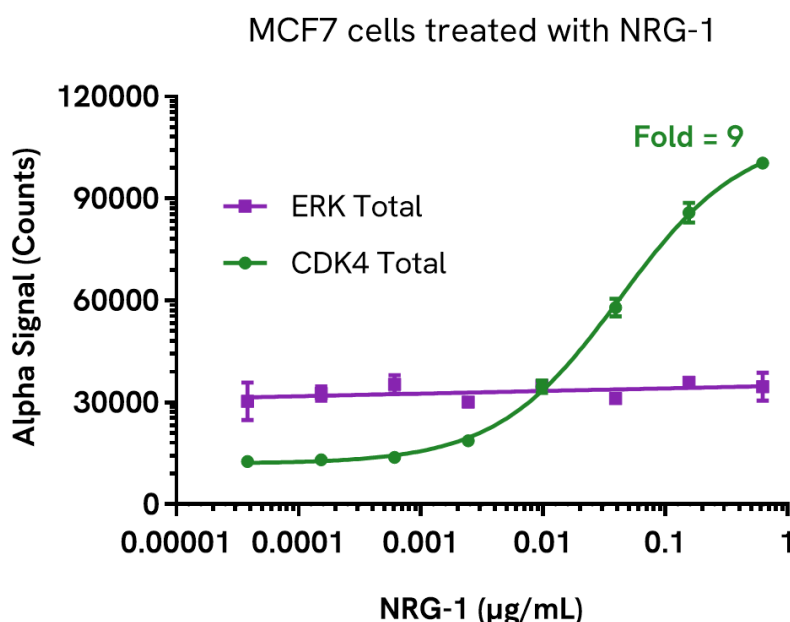

Data obtained with a 2-plate, 2-incubation protocol. RT4 cells were seeded at 40K cells/well in a 96-well plate and incubated overnight. Cells were treated with increasing concentrations of BSJ-03-204 (CDK4/6 PROTAC) for 24 hours. Cells were lysed with Lysis Buffer and assayed separately for Total CDK4 and ERK using respective SureFire Biotin Free kits. Equivalent to approximately 4,000 cells/datapoint.

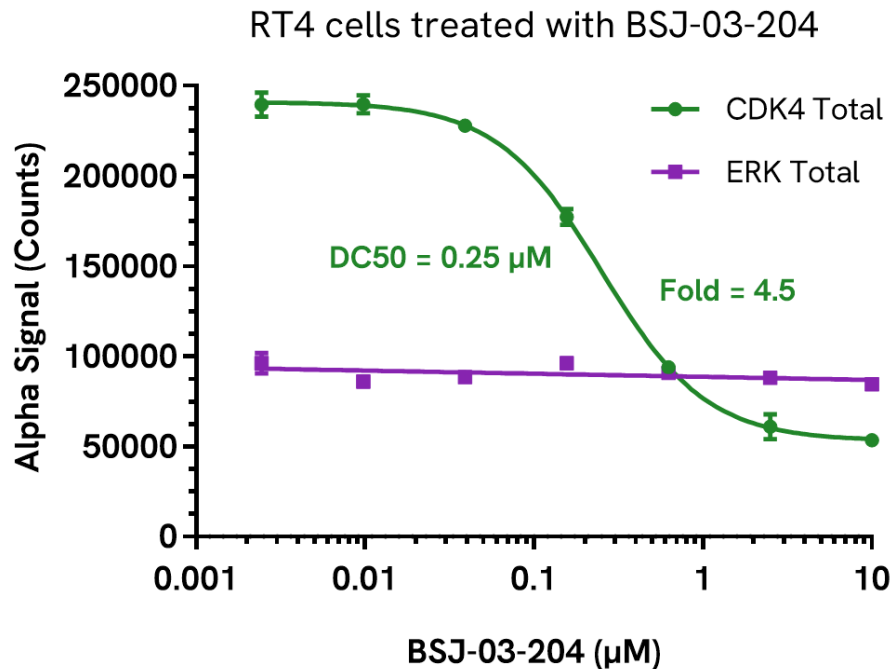

The information provided in this document is for reference purposes only and may not be all-inclusive. Revvity, Inc., its subsidiaries, and/or affiliates (collectively, "Revvity") do not assume liability for the accuracy or completeness of the information contained herein. Users should exercise caution when handling materials as they may present unknown hazards. Revvity shall not be liable for any damages or losses resulting from handling or contact with the product, as Revvity cannot control actual methods, volumes, or conditions of use. Users are responsible for ensuring the product's suitability for their specific application. REVVITY EXPRESSLY DISCLAIMS ALL WARRANTIES, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, REGARDLESS OF WHETHER ORAL OR WRITTEN, EXPRESS OR IMPLIED, ALLEGEDLY ARISING FROM ANY USAGE OF ANY TRADE OR ANY COURSE OF DEALING, IN CONNECTION WITH THE USE OF INFORMATION CONTAINED HEREIN OR THE PRODUCT ITSELF"

