

AlphaLISA™ SureFire® Biotin Free

Human and Mouse Histone H3 Total Detection Kit

Product number: ASBF-THISH3-A500, ASBF-THISH3-A10K,
ASBF-THISH3-A50K, ASBF-THIS3-A-HV



Kit specificity:


This assay kit contains antibodies which recognize distinct epitopes on Histone H3 protein. The protein detected by this kit corresponds to UniProt ID P68431. Histone H3 is also known as Histone H3.1. These antibodies recognize Histone H3 of human and mouse origin. Other species should be tested on a case-by-case basis.

Control lysate information:


Positive Control Lysate: Prepared from U2OS cells cultured to confluence in T175 flasks in 10% FBS containing medium, then treated with 10 μ M Nocodazole for 16 hours, and lysed with 120 mL of Lysis Buffer C. The Histone H3 Total assay is compatible with Lysis Buffer B.

Representative data:


Data obtained with a 2-plate, 2-incubation protocol. HeLa cells were seeded at 100K cells/well in a 96-well plate and treated with increasing concentrations of Trichostatin A (TSA) for 4 hours. Cells were lysed with 200 μ L Lysis Buffer B for 30 minutes with shaking. Lysate was further diluted with Lysis Buffer B and assayed separately for Acetyl Histone H3 (K27) and Histone H3 Total using respective SureFire Biotin Free kits. Equivalent to approximately 100 cells/datapoint.

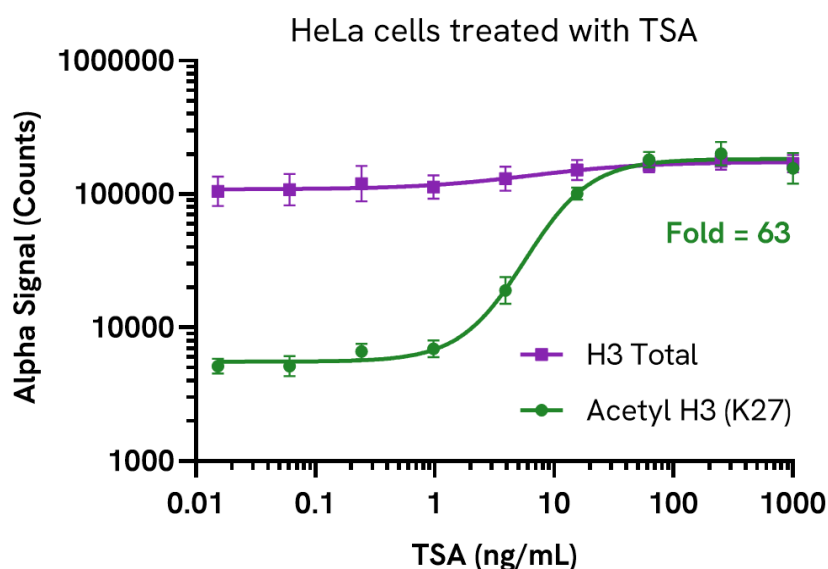


Data obtained with a 2-plate, 2-incubation protocol. Karpas 299 cells were seeded at 400K cells/well in a 96-well plate and treated with increasing concentrations of A485 for 4 hours. Cells were lysed with the addition of 50 μ L of 5X Lysis Buffer Lysis Buffer B for 30 minutes with shaking. Lysate was assayed separately for Acetyl Histone H3 (K27) and Histone H3 Total using respective *SureFire* Biotin Free kits. Equivalent to approximately 16,000 cells/datapoint.

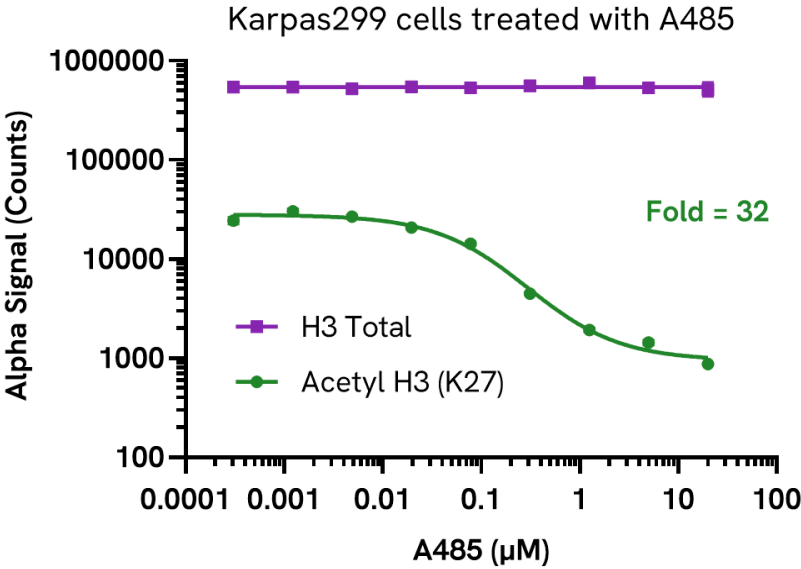

Data obtained from measurement of Histone H3 Total in various cell types and lysed with Lysis Buffer B for 30 minutes with shaking. Equivalent to approximately 500 cells/datapoint for adherent cells or 1,000 cells/datapoint for suspension cells.

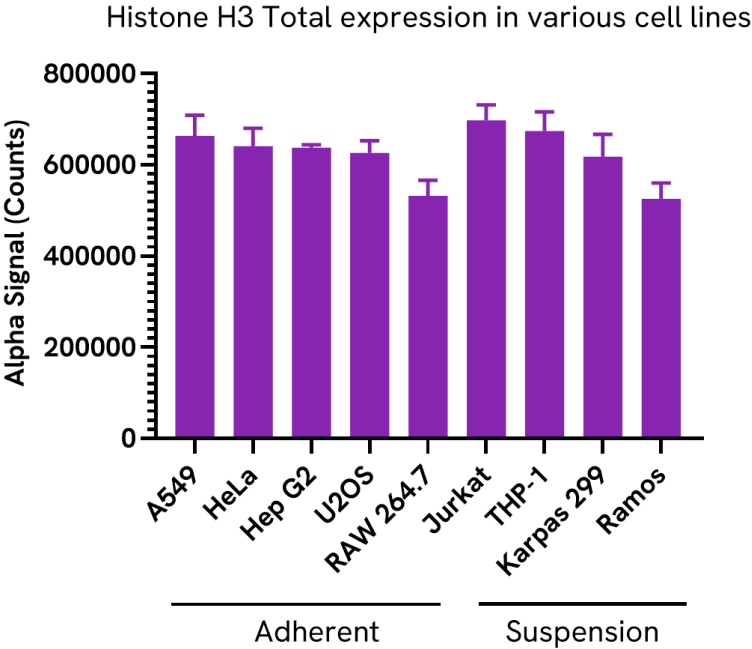

Data obtained from measurement of Histone H3 Total in PBMCs lysed with Lysis Buffer B for 30 minutes with shaking. Lysate was serially diluted with Lysis Buffer B and assayed using the Histone H3 Total *SureFire* Biotin Free kit.

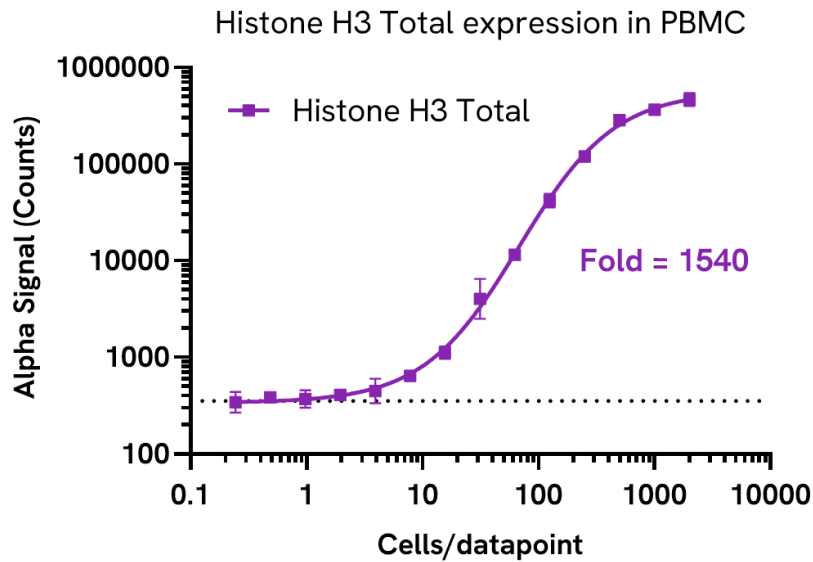

The information provided in this document is for reference purposes only and may not be all-inclusive. Revvity, Inc., its subsidiaries, and/or affiliates (collectively, “Revvity”) do not assume liability for the accuracy or completeness of the information contained herein. Users should exercise caution when handling materials as they may present unknown hazards. Revvity shall not be liable for any damages or losses resulting from handling or contact with the product, as Revvity cannot control actual methods, volumes, or conditions of use. Users are responsible for ensuring the product's suitability for their specific application. REVVITY EXPRESSLY DISCLAIMS ALL WARRANTIES, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, REGARDLESS OF WHETHER ORAL OR WRITTEN, EXPRESS OR IMPLIED, ALLEGEDLY ARISING FROM ANY USAGE OF ANY TRADE OR ANY COURSE OF DEALING, IN CONNECTION WITH THE USE OF INFORMATION CONTAINED HEREIN OR THE PRODUCT ITSELF”

